

Trucuff LN Larraged Washington

Truc Cuff TI Tracked Tube Intelion SINCLEUSE

NBSK HUMANON SINGLEUSE

Tru-Cuff MLM&TT



The Tru-Cuff[™] indicates pressure in endotracheal tube and laryngeal mask cuffs, minimizing the potential for airway trauma from over inflation of the cuff. The Tru-Cuff[™] allows the clinician to easily inflate airways to a clinically recommended pressure with a single patient device. Easy to read color indicator bands provide you with at-aglance feedback indicating cuff pressure, therefore improving patient safety and minimizing the risk of airway damage. The Tru-Cuff[™] is available with the complete family of AES Ultra laryngeal masks or can be purchased as a stand-alone for laryngeal masks and endotracheal tubes.



AES, Inc.


Specifications and ordering information

Tru-Cuff [™] LM + TT Ordering Info				
Size	Item #	Description	Packaging	
LM	105004060	Laryngeal Inflation Syringe	Box of 50	
TT	105002030	Endotracheal Inflation Syringe	Box of 50	
Com	ponent	Material		

Syringe Body	Polycarbonate
Bellows	Silicone
Plunger	Isoprene

Presure Range

Tru-Cuff [™] LM	40 – 60 cmH ² O
Tru-Cuff [™] TT	20 – 30 cmH ² O

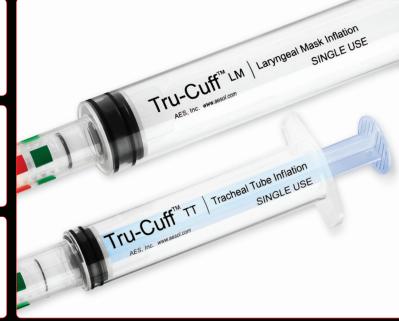

Distributed by:



Call toll free: 800.426.5007 or visit: www.aesol.com AES, Inc.



Technical data may be modified without further notice.

AES, Inc • 24301 Roberts Drive • Black Diamond, WA 98010